


## Fluorescence Microscopy Filter Cube Schematic Diagram

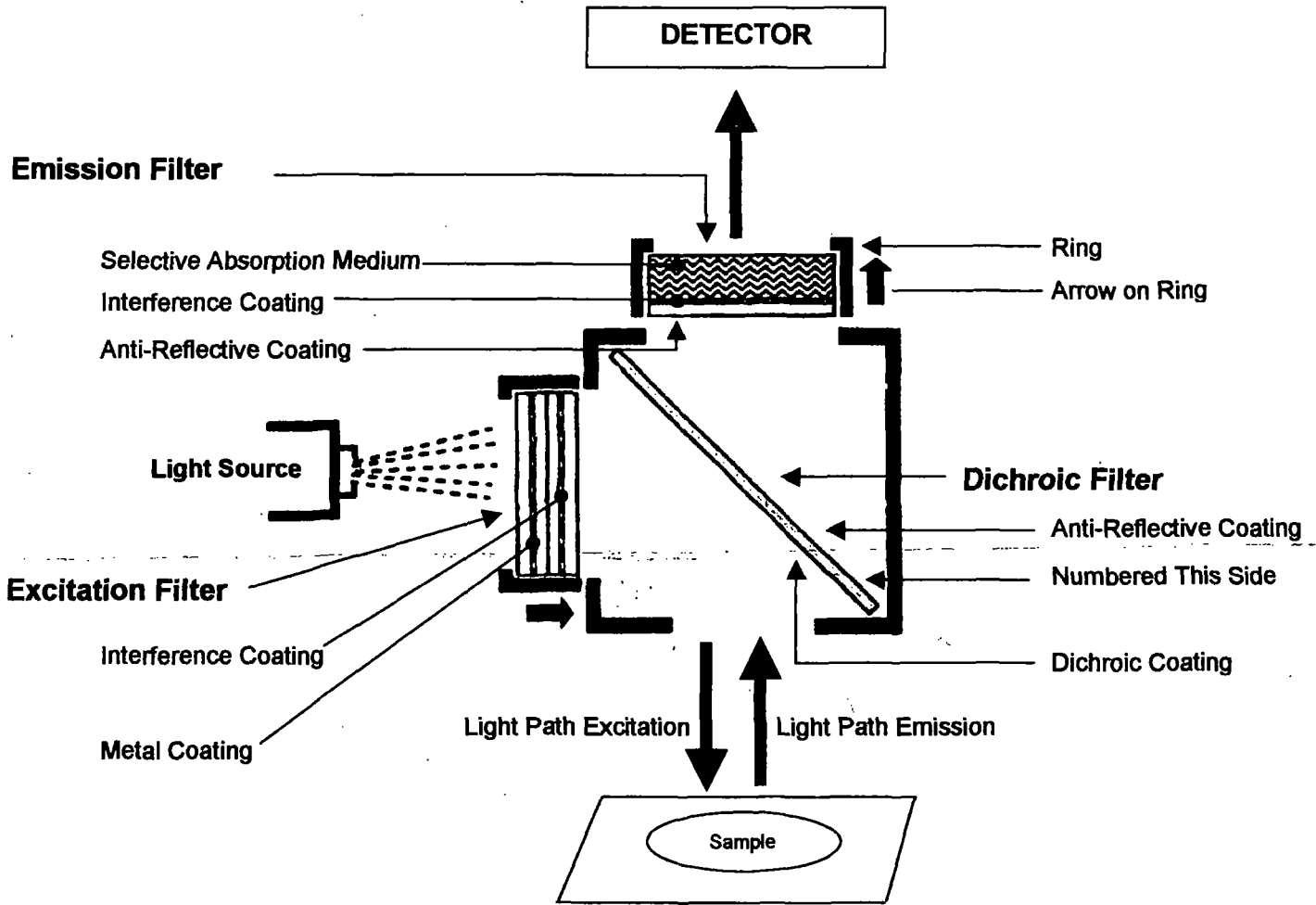

### Important Notes:

#### Orientation of Filters:

Omega exciter & emission filters have an arrow on ring pointing in the direction of the light path.

#### Care & Handling of Filters:

Precision optical components should be handled with reasonable care. Hold by edges only.

#### Cleaning of Optical Components:

Remove any foreign particles with a puff of dry air. Wipe gently with a soft, lint-free cloth. A final wipe with a few drops of pure anhydrous alcohol will result in a clean, undamaged component.